



 borer medical

The NEW Borer PCDS Premium -
because quality matters!

 **borer**
advanced cleaning solutions



Innovation and quality from Borer

A high quality and innovative system. This was the goal for the development of our new conveyor system for process chemicals. On the one hand, quality is more important than ever in Central Sterile Supply Departments (CSSD) and the instrument reprocessing cycle¹. On the other, today's customers demand sophisticated yet simple solutions.

Main functions

- Continuous, uninterrupted supply of process chemicals to the connected WDs (washer-disinfectors).
- Simpler, better logistics aimed at improving the efficiency of the cleaning and disinfecting process.

Main advantages

- Constant supply to WDs
The connected WDs are regularly supplied with process chemicals without suspension. This avoids the unnecessary suspension of a washing cycle just because a container is empty.
- Uses large containers
The PCDS can hold large 220 litre containers for example. This eliminates the need for changing smaller containers, which is laborious and time-consuming.

- Supply to several WDs at once
The new PCDS can supply up to 10 WDs with process chemicals at the same time, regardless of their location.
- Contactless replacement of containers
The lifting mechanism means that staff do not have to touch the suction lances, which protects them on the one hand and makes handling easier on the other.

Unique properties

- Notification to change the container
The system sends email notification as soon as the container reaches the preset (critical) fill level. The system monitors itself so you do not have to.
- Integrated web server
The system and the system functions are easily monitored from a remote PC.

Simple to use and time-saving

The new Borer PCDS Premium meets the needs of CSSD staff: the user-friendly system reduces handling and thus helps save time and money. Furthermore, it guarantees absolute reliability in the initial stages of a successful cleaning and disinfection process.


Just perfect for CSSD Managers

As a CSSD manager you are responsible for a strict, efficient and smooth cleaning and disinfection process.

You also ensure the protection of staff and patients with the dedicated products and services. That is why quality is extremely important and why the new Borer PCDS Premium is just perfect. Because quality matters!

We would be pleased to present the new Borer PCDS Premium to you in more detail. For more information please contact our distributor partner in your area or visit www.borer.ch/en/borer-pcds-premium.

¹ Fig. The instrument reprocessing cycle

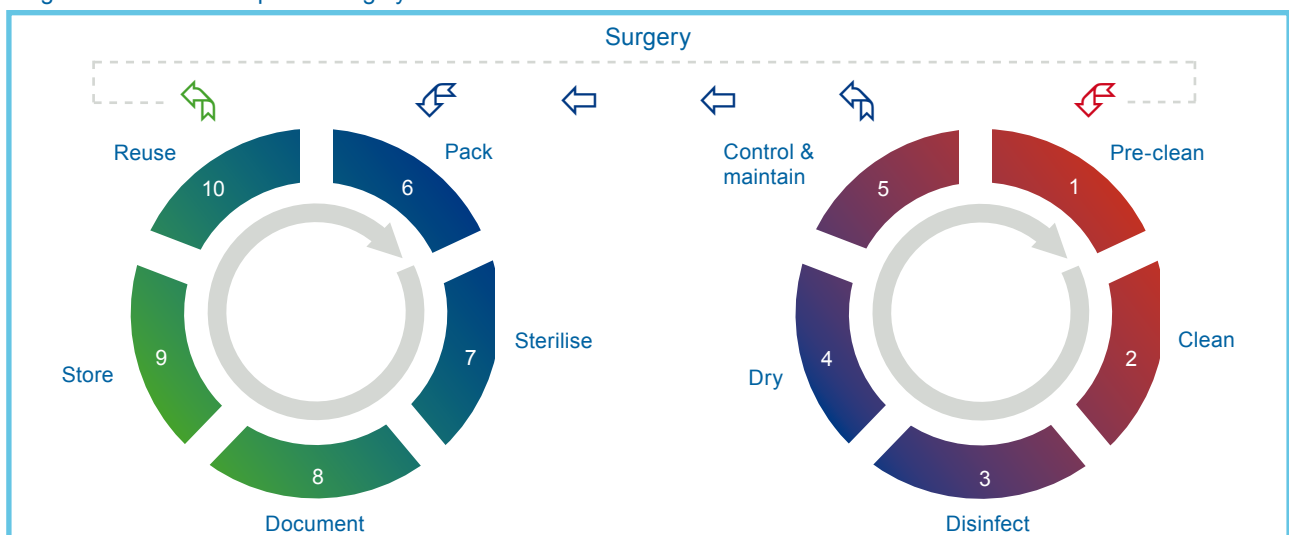

Borer Chemie AG

Gewerbstrasse 13, 4528 Zuchwil / Switzerland
Tel. +41 32 686 56 00, Fax +41 32 686 56 90
www.borer.ch, office@borer.ch

deconex[®]